
REMOTE | **Getting Started**



Usage: The MysteryVibe Remote allows you to change intensity and vibrations on any of your MysteryVibe devices from a distance of up to 3 metres (10 feet), and can control up to 3 MysteryVibe devices to stimulate multiple areas at the same time.

Connecting: Turn on Bluetooth on the devices you want to use with the Remote. Press and hold the Remote's Scan button. Each channel light will flash yellow when connected to a device, and turn green to indicate which channel is currently under control.

Charging: Connect the USB wire to the port on the bottom of your Remote and the other end to any USB plug for power. The status light will blink green to indicate charging and turn green when charging is complete. Charge for 60 minutes before use.

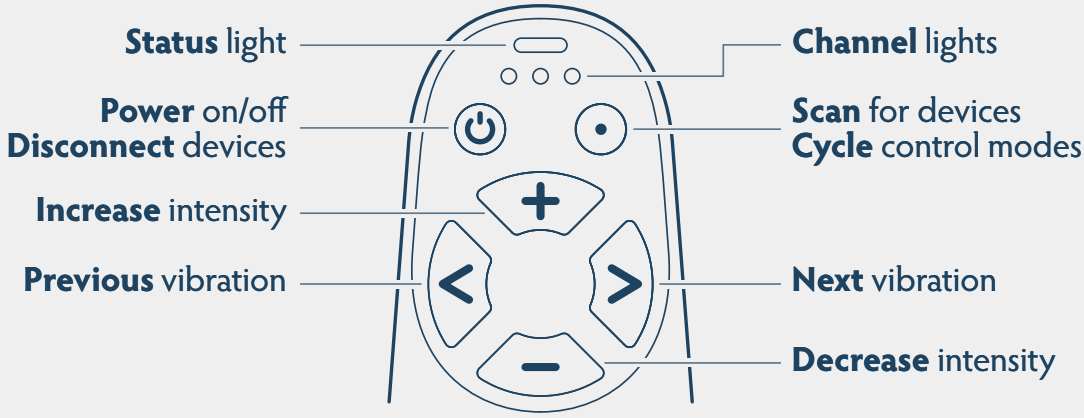

Warranty information: for registration help, email support@mv.health. Your seller can confirm the process and terms for your purchase. B2B support copy.

FDA-registered Class II therapeutic genital vibrator accessory (21 CFR § 884.5960, Product Code KXQ). Not intended to diagnose, treat, cure, or prevent any disease or condition.

6905-2007-100-1.6