

LEGATO | **Getting Started**

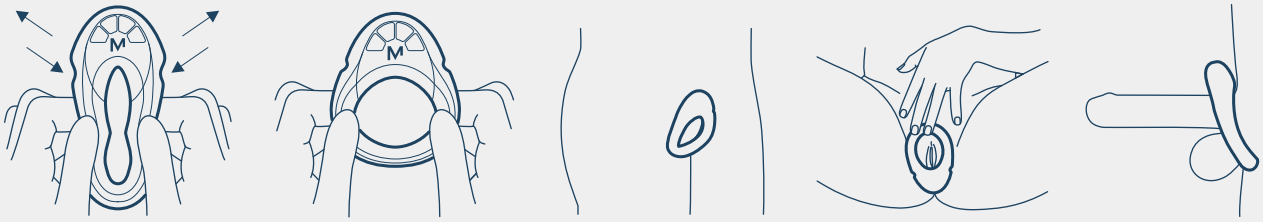


Usage: Turn Legato on and stretch/squeeze the sides inwards/outwards to take the shape that works best for you. Hold Legato against your body & feel its full 360° vibrations.

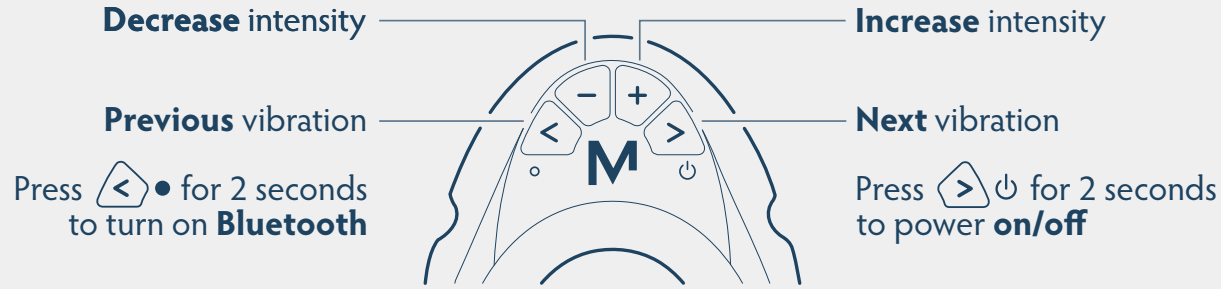

Charging: Connect the USB wire to the port on the top of your Legato and the other end to any USB plug for power. The light on Legato will blink amber to indicate charging. When the light turns green, charging is complete. Charge for 60 minutes before use.

App: Download the free MysteryVibe App from the Apple App Store or Google Play. Connect and use the App as a remote control for your Legato. Download different vibration patterns or create your own.

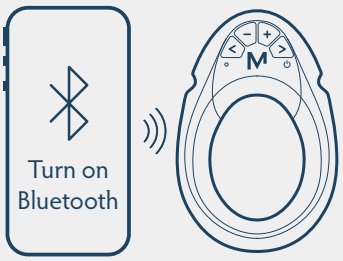

FDA-registered Class II therapeutic genital vibrator (21 CFR § 884.5960, Product Code KXQ). Not intended to diagnose, treat, cure, or prevent any disease or condition

6914-2201-122-1.2