

MOLTO | **Getting Started**




Usage: Molto is designed to mimic the width of one finger, adapt to your anus or vagina and stimulate precisely where needed. Gently and slowly insert Molto to reach and stimulate the anus, prostate or vagina.

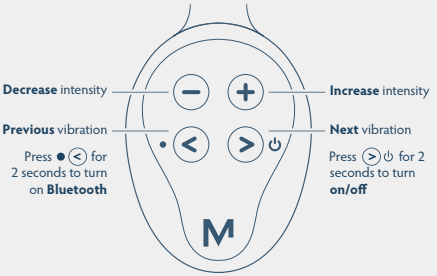

Charging: Connect the USB wire to the port on the bottom of your Molto and the other end to any USB plug for power. The light on Molto will blink amber to indicate charging. When the light turns green, charging is complete. Charge for 60 minutes before use.

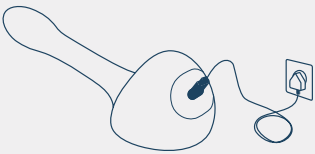

App: Download the free MysteryVibe App from the Apple App Store or Google Play. Connect and use the App as a remote control for your Molto. Download different vibration patterns or create your own.

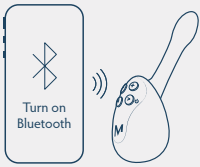

Warranty information: for registration help, email support@mv.health. Your seller can confirm the process and terms for your purchase. B2B support copy.

6915-2207106-14

FDA-registered Class II therapeutic genital vibrator (21 CFR § 884.5960, Product Code KXQ). Not intended to diagnose, treat, cure, or prevent any disease or condition