

CRESCENDO 2 | **Getting Started**




Usage: Crescendo 2 is designed to mimic two fingers so you can easily bend & adapt it to take the shape that works best for you and deliver targeted vibrations, precisely where you need it.

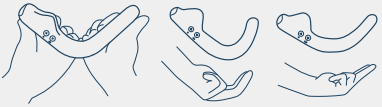

Charging: Connect the USB wire to the port on the back of your Crescendo 2 and the other end to any USB plug for power. The light on Crescendo 2 will blink amber to indicate charging. When the light turns green, charging is complete. Charge for 60 minutes before use.

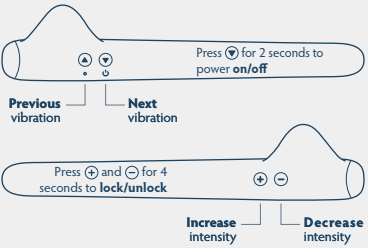

App: Download the free MysteryVibe App from the Apple App Store or Google Play. Connect and use the App as a remote control for your Crescendo 2. Download different vibration patterns or create your own.

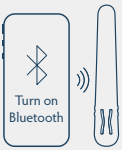

Warranty information: for registration help, email support@mv.health. Your seller can confirm the process and terms for your purchase. B2B support copy.

FDA-registered Class II therapeutic genital vibrator (21 CFR § 884.5960, Product Code KXQ). Not intended to diagnose, treat, cure, or prevent any disease or condition

6908-2103-107-13