


# TENUTO 2 | **Getting Started**

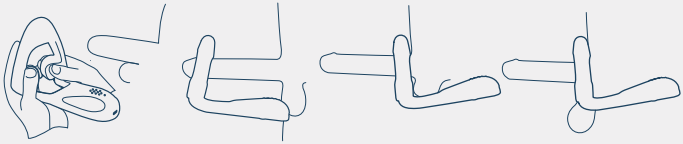


**Usage:** Apply lube to both Tenuto 2 and your penis. Turn Tenuto 2 on, twist & stretch the front outwards and place it on the head of the penis. Once aroused, slide it down the penis shaft so it sits at the base. Experience powerful vibrations on the penis, perineum and stimulate your partner during intercourse.

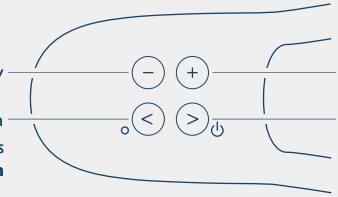


**Charging:** Connect the USB wire to the port on the back of your Tenuto 2 and the other end to any USB plug for power. The light on Tenuto 2 will blink amber to indicate charging. When the light turns green, charging is complete. Charge for 60 minutes before use.

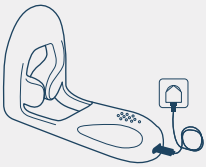

**App:** Download the free MysteryVibe App from the Apple App Store or Google Play. Connect and use the App as a remote control for your Tenuto 2. Download different vibration patterns or create your own.

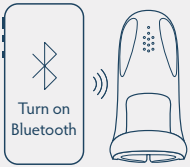


Warranty information: for registration help, email support@mv.health. Your seller can confirm the process and terms for your purchase. B2B support copy.

FDA-registered Class II therapeutic genital vibrator (21 CFR § 884.5960, Product Code KXQ). Not intended to diagnose, treat, cure, or prevent any disease or condition

6909-2103-106-13